

panakeia®



Two Gas System
Basic



Three Gas System
Advanced



Panakeia Medtech Pvt Ltd

Factory Address : AMTZ Campus, VM Steel Project,
Pragathi Maidan, Visakhapatnam - 530 031, (AP)
Mob.: 9811340469, 9811340468

panakeia®

PANAKEIA 'PISCES' ANAESTHESIA WORKSTATION

→ **ANAESTHESIA WORKSTATION** is equipped with integrated Anaesthesia Ventilator & Provision for two Thermo/Pressure Compensated Flow independent selectatec mount Vaporizer.

→ Structure

Made of Aluminum Extruded Tubes with Powder Coated Finish.

Drawers fitted with sliding channels for smooth movement.

Sliding writing tray at middle table for use.

Four Nos. big diameter Castor Wheels. Front Castors are equipped with individual brakes.

Handle fitted in both the side for easy maneuverability.

Space for mounting up to two nos. vaporizers- Halothane, Isoflurane, Sevoflurane.


→ Gas inlet

Cylinder Inlet: Four nos., Pin indexed Yoke, each for O₂ & N₂O. Yokes are identified by clear marking.

Central Pipe line Inlet: Non- interchangeable pipeline connection for O₂, N₂O, Air



Advanced

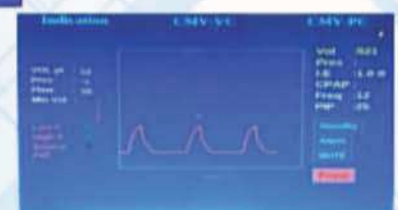

Basic

→ **Pressure Regulators: 2 stage regulation**

Each gas service is fitted with high performance Regulators for two stage reduction.

→ **Pressure Indicators (Pressure Gauges): High & Low Pressure**

The Pressure gauge Panel is mounted with big pressure gauges for individual cylinder pressure indication of O₂ & N₂O. The gauges are identified by the gas service marked on them as well as by Colour coding.

→ **Flowmeter with Oxygen Ratio Controller:**

It is equipped with a Cascade type Five Tube Rotameter Assembly for flow control of O₂, N₂O & Air. It is fitted with high and low range tubes each of O₂ & N₂O for precise control of gas flow.

GasFlowmeterRange O₂Double –

0.1 1 litre/min 1-10 litre/min N₂O Double 0.1 – 1 litres/min

Gas	Flowmeter	Range
O ₂	Double	0.1 – 1 litre/min 1-10 litre/min
N ₂ O	Double	0.1 – 1 litres/min 1 – 12 litres/min
Air	Single	200 – 15000 cc/min



→ **Oxygen Low Pressure Alarm:**

This safety device generates an audio & visual alarm when oxygen pressure in the rotameter falls below 500-600ml.

→ **Air/N₂O Interlock:**

This interlock system eliminates the possibility of selecting Air with N₂O.

→ **Vaporizer Manifold:**

It is provided with Two Thermo-Compensated Vaporizer of user's choice. Space provided for mounting another Vaporizer,

→ **Oxygen Flush:**

This flush unit is designed to supply emergency oxygen to the patient at the rate of 30 – 75 lpm. The push button switch of the unit is mono stable type and to be pressed for use.

→ Anaesthesia Ventilator

Micro-controller based Electronically controlled Pneumatic driven Ascending Bellow Ventilator

7 Inches Touch Screen

Compatible with Circle Absorber, Bain's circuit

Comprehensive alarms {Patient disconnection & High pressure, Low Insp. Pressure}

Modes : CMV (VC), CMV (PC)

Modes (Optional) : SIMV (VC), SIMV (PC) WITH PSV

Tidal Volume : 50 to 1500 ml.

Frequency : 6 to 40 bpm

Inspiratory time : 0.4 to 3 seconds

I/E ratio (variable) : 4:1 to 1:4

Inspiratory hold : 0 to 2 seconds

PIP : 10 to 60 cm of water

PEEP/CPAP : 0 to 20 cm of water



→ Standard Accessories:-

Circle Absorber System

One Thermo-Compensated Vaporizer

Bain's Circuit

Pediatric Anaesthesia Circuit

Colour Coded High Pressure Tubing of O₂, N₂O, Air

Face Mask of 2 Sizes ;

2 Nos. Rebreathing Bags;



CE ISO 9001:2015
ISO 13485:2016

